



# FenomPro<sup>®</sup>

Asthma Monitor

## Instructions for Use



**R<sub>x</sub> Only**

# Table of Contents

<b>Definitions</b>	<b>4</b>		
<b>Chapter 1: System Overview</b>	<b>5</b>		
System Description	5		
Indications for Use	5		
Clinical Limitations	5		
Risks to Health	5		
Fenom Pro Components	6		
Table 1: Fenom Pro Components	6		
Display Buttons	7		
Table 2: Button and Indication Icons	7		
<b>Chapter 2: Safety and Warnings</b>	<b>8</b>		
Safety Instructions	8		
Compliance	8		
Warnings	8		
Electromagnetic Emissions	9		
Electromagnetic Immunity	9		
<b>Chapter 3: Fenom Pro Quick Start Guide</b>	<b>10</b>		
<b>Chapter 4: FeNO Measurement Preparation</b>	<b>12</b>		
Wake up Device	12		
Pre-test Check	12		
<b>Chapter 5: Perform FeNO Measurement</b>	<b>13</b>		
Perform a FeNO Test	13		
Remove Mouthpiece	14		
<b>Chapter 6: Quality Control</b>	<b>15</b>		
Perform a QC Test	15		
QC User Status	15		
Create QC Users	15		
<b>Chapter 7: Power Off Device</b>	<b>16</b>		
<b>Chapter 8: Practice Mode</b>	<b>16</b>		
<b>Chapter 9: Device Setup</b>	<b>17</b>		
Initial Setup	17		
Configuration Settings	17		
Device Settings	17		
Time/Date	17		
Language	17		
		System Information	19
		Test Incentive Sound	19
		Add Tests	19
		Quality Control Status	19
		<b>Chapter 10: General Care</b>	<b>20</b>
		Operating Conditions	20
		Cleaning	20
		Handling	20
		Storage	20
		Preventive Inspections	21
		Rechargeable Battery	21
		Maintenance	21
		Disposal of Used/Expired Equipment and Consumables	21
		Limited Warranty	21
		<b>Chapter 11: Troubleshooting</b>	<b>22</b>
		Support	22
		Add Tests	22
		Error and Codes	23
		Table 4: Error Codes	23
		<b>Chapter 12: Technical Data</b>	<b>26</b>
		Dimensions and Weight	26
		Electrical Data	26
		Exhaled NO Performance	26
		Linearity	26
		Precision	26
		Accuracy	26
		Limit of Detection	26
		Measurement Range	26
		Exhalation Parameters	26
		<b>Chapter 13: Device Performance</b>	<b>27</b>
		Accuracy	28
		Limit of Detection	28
		<b>Chapter 14: Reference</b>	<b>29</b>
		Symbol Explanation	29
		<b>Chapter 14: Parts and Accessories</b>	<b>30</b>
		Parts	30
		Accessories	30
		<b>Bibliography</b>	<b>31</b>

## Definitions

The following are the definitions for terms and abbreviations used in this manual

<b>FeNO</b>	Fractional exhaled nitric oxide – Amount of nitric oxide in the exhaled breath originating from the bronchial passages, not the nasal passages nor upper airway.
<b>FEV1</b>	Forced Expiratory Volume in One Second - Volume of air that can be forcibly exhaled from the lungs in the first second of a forced expiratory maneuver, measured in liters.
<b>FEV6</b>	Forced Expiratory Volume in Six Seconds - Volume of air that can be forcibly exhaled from the lungs in the six seconds of a forced expiratory maneuver, measured in liters.
<b>FVC</b>	Forced Vital Capacity - After the patient has taken in the deepest possible breath, this is the volume of air that can be forcibly and maximally exhaled out of the lungs until no more can be expired, usually measured in liters.
<b>NO</b>	Nitric oxide – Produced by the human lung and is present in the exhaled breath. It has been implicated in the pathophysiology of lung diseases, including asthma.
<b>PEF</b>	Peak Expiratory Flow - Maximal flow (or speed) achieved during the maximally forced expiration initiated at full inspiration, measured in liters per minute or in liters per second.
<b>Spirometry</b>	Common office test used to assess how well a patient's lungs work by measuring how much air is inhaled, how much is exhaled, and how quickly it is exhaled.

## Chapter 1: System Overview

### System Description

SpiroSure Fenom Pro™ is a point-of-care breath analyzer that uses electrochemical sensor technology to measure the fraction of exhaled nitric oxide (FeNO), a marker for airway inflammation, in human exhaled breath. Measurement of FeNO by Fenom Pro is quantitative, non-invasive (only the mouthpiece comes in contact with the patient), simple and safe. Fenom Pro is designed as a hand-held device for measuring FeNO in exhaled breath from humans. The level of exhaled nitric oxide (NO) is frequently increased in some inflammatory processes such as asthma. The fractional NO concentration in expired breath can be measured by the Fenom Pro™ device according to guidelines for NO measurement established by the American Thoracic Society (ATS) and European Respiratory Society (ERS) [1].

Fenom Pro™ provides direct sampling with result report within 30 seconds of sequentially collected and analyzed exhaled air. No subsequent specific specimen collection, specimen preparation, or reagents are required. The emissions characteristics of the Fenom Pro™ device make it suitable for use in point of care sites and hospitals (CISPR 11 class A).

### Indications for Use

Fenom Pro Nitric Oxide Test is a portable, non-invasive device to measure fractional exhaled nitric oxide (FeNO) in human breath. FeNO is increased in some airway inflammatory processes, such as asthma, and often decreases in response to anti-inflammatory treatment. Measurement of FeNO by Fenom Pro is a method to measure the decrease in FeNO concentration in asthma patients that often occurs after treatment with anti-inflammatory pharmacological therapy as an indication of therapeutic effect in patients with elevated FeNO levels. FeNO measurements are to be used as an adjunct to established clinical assessments. Fenom Pro is suitable for children, approximately 7-17 years, and adults 18 years and older.

Testing using the Fenom Pro should only be done in a point-of-care healthcare setting under professional supervision. Fenom Pro should not be used in critical care, emergency care or in anesthesiology.

### Clinical Limitations

Fenom Pro may not be used by children under the age of approximately 7 years, including infants, as measurement requires patient cooperation. Fenom Pro may not be used by children under the age of 7 years, or by patients who are unable to understand and execute the instructions given by healthcare providers, as measurement requires patient cooperation.

Fenom Pro should not be used in critical care, emergency care, or in anesthesiology.

All subjects should refrain from eating or drinking for at least 60 minutes before the FENOM test. Recent intake of nitrate rich food, such as lettuce, spinach, beets, walnuts, peanuts, and animal organs, can lead to increased FeNO levels [2].


Smoking reduces exhaled NO levels [5]. Fenom results obtained from subjects who smoke should only be considered after considering the subject's smoking history and the potential impact on NO levels.

### Risks to Health

There are no known direct risks to patient health posed by use of Fenom Pro. However, failure of the test to perform as indicated or erroneous interpretation of results may lead to improper patient management.

Therefore, use of FeNO measurement results to adjust a treatment regimen without consideration of other clinical factors could pose a risk.

## Fenom Pro Components

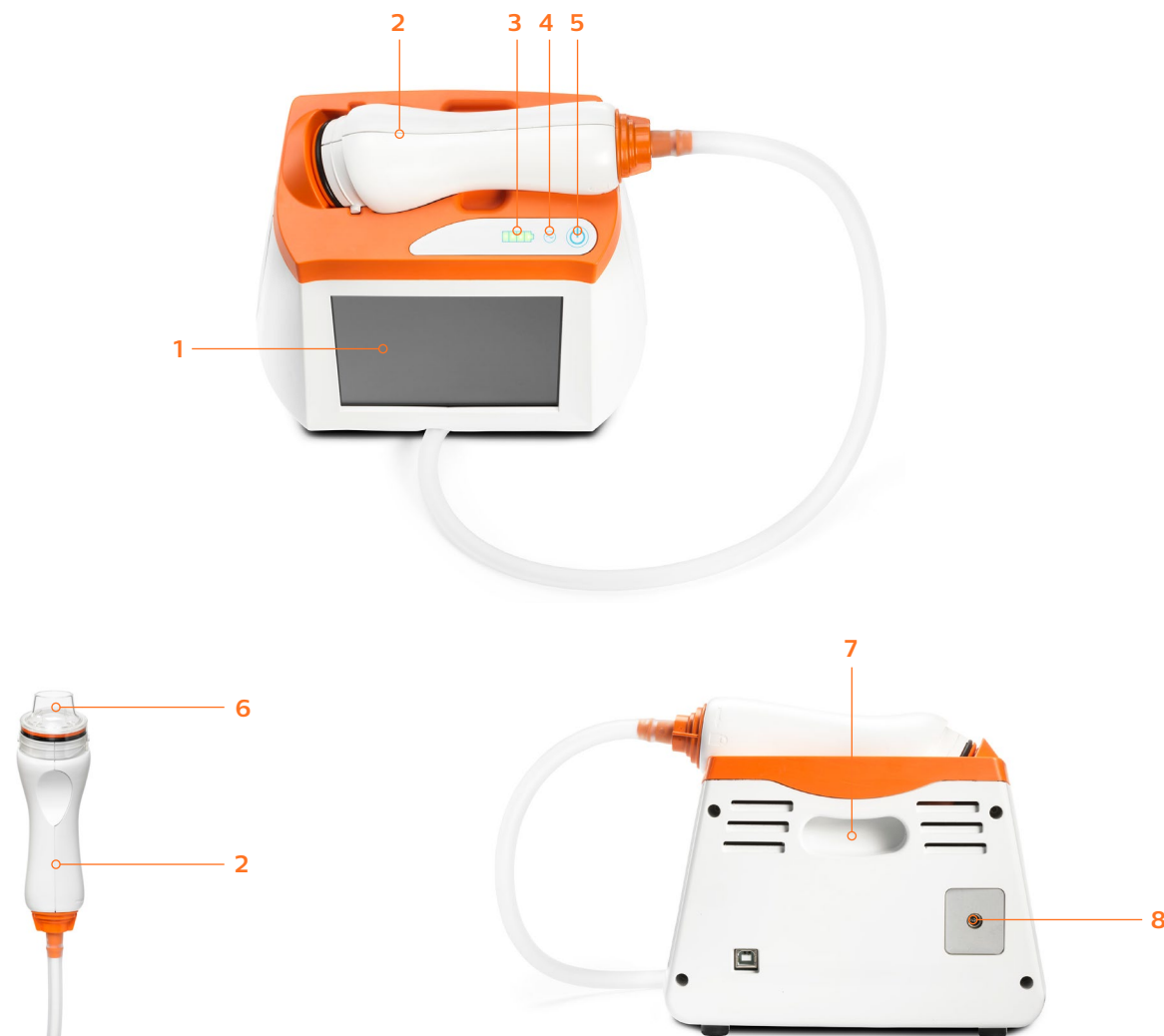

**Table 1: Fenom Pro Components**

No.	Description
1	Touch Screen
2	Handpiece
3	Battery Indicator –Battery strength is below 25% if only one bar is illuminated
4	AC Power Indicator - Indicator is <b>green</b> when the device is powered on and connected to an electrical outlet.
5	Power Button – Hold for one second to power on/off.
6	Single-Patient-Use Mouthpiece (accessory)
7	Carrying Handle
8	24 V Power Connection

## Display Buttons

There are several button icons that Fenom Pro utilizes to help you easily navigate through the menu screens.

**Table 2: Button and Indication Icons**

Buttons Icon	Name	Description
	<b>Settings Button</b>	Button used to open the Settings Menu. This menu allows for setting Time/Date, selecting Language, viewing System Information, selecting Volume Level, and Ordering Tests.
	<b>Test License Status Button (Tests available)</b>	Button used to open Order Tests box. <b>Green</b> checkmark indicates that the device has tests available.
	<b>Test License Status Button (Few tests available)</b>	Button used to open Order Tests box. <b>Red</b> exclamation point indicates few tests are remaining. Contact your distributor to order additional tests.
	<b>Quality Control Status Button (QC Status - Pass)</b>	Button used to indicate current QC Mode Status. <b>Green</b> icon indicates current QC Mode status is passed.
	<b>Quality Control Status Button (QC Status - Fail)</b>	Button used to indicate current QC Mode Status. <b>Red</b> icon indicates current QC Mode status is failed or expired.
	<b>Quality Control Status Button (QC Status - Disabled)</b>	Button used to indicate current QC Mode Status. Gray icon indicates current QC Mode status is temporarily disabled or analyzing.

## Chapter 2: Safety and Warnings

### Safety Instructions

The following safety instructions apply in the handling and operation of Fenom Pro:

- DO NOT inhale through the device.
- DO NOT inhale through the mouthpiece.
- DO NOT exhale beyond the limits of your physical ability.
- Discontinue measurements if the breath maneuver is laborious for the patient.
- DO NOT allow use of Fenom Pro within 15 minutes after performing spirometry testing such as: FEV1, FEV6, FVC, PEF, etc.
- DO NOT allow use of Fenom Pro within 60 minutes after exercising or smoking.
- DO NOT use the Fenom Pro device without a new single-patient-use mouthpiece.
- DO NOT perform more than six breath attempts on a single patient within one day.
- If repeated measurements taken for at least two consecutive days are more than 10 ppb different contact Spirosure, Inc. Technical Support.
- DO NOT bring Fenom Pro in a room containing magnetic resonance equipment.
- DO NOT bring Fenom Pro to a room adjacent to magnetic resonance equipment.

### Compliance

Fenom Pro is CE-marked according to the *In Vitro Diagnostics Directive 98/79/EC*.

Fenom Pro is RoHS compliant according to *Directive 2011/65/EU Restriction of Hazardous Substances in Electrical and Electronic Equipment*.

### Warnings

The following warnings apply in the handling and operation of Fenom Pro:

- Fenom Pro should only be operated by trained healthcare professionals.
- Operate Fenom Pro as stated in this manual. Spirosure accepts no responsibility for damaged equipment or faulty results if the equipment is not handled according to this manual.
- DO NOT use a damaged Fenom Pro device, damaged components, or damaged accessories.
- Use only power supply unit provided.
- Keep the device out of water. Ensure no liquid is spilled or dripped on the device.
- DO NOT use the Fenom Pro device adjacent to or stacked with other equipment because it could result in improper operation.

- DO NOT block device vents and ports while in use or while charging.
- DO NOT drop the device or subject it to strong impact.
- No modification of the Fenom Pro device, handpiece, or mouthpiece is allowed.
- DO NOT use Fenom Pro in the proximity of areas where volatile substances such as organic fluids or disinfectants are being used. Special attention should be paid to aerosols and disinfection baths.
- DO NOT use Fenom Pro in the presence of flammable vapors or liquids
- Use of substances containing alcohol close to the Fenom Pro device may cause erroneous measurement results. [3]
- Single-patient-use mouthpiece should be used immediately after opening.
- DO NOT reuse the single-patient-use mouthpiece on other patients.
- When not in use, the Fenom Pro device should be stored in the packaging provided. (See Chapter 12: General Care.)
- DO NOT open, crush, heat above 140 °F/60 °C, or incinerate the lithium-ion battery in the device.

### Electromagnetic Emissions

The emissions characteristics of this equipment make it suitable for use in hospitals and other healthcare settings (CISPR 11 class A).

### Electromagnetic Immunity

Fenom Pro has been tested to comply with the emission and immunity requirements described in IEC 60601-1-2:2014 (4th Edition) General requirements for basic safety and essential performance - Collateral standard: Electromagnetic compatibility - Requirements and tests and AIM 7351731:2017 Medical Electrical Equipment and System Electromagnetic Immunity Test for RFID Readers.

Fenom Pro is not compatible with magnetic resonance systems, and is labeled as MR Unsafe. The Fenom Pro should not be used in a room containing a magnetic resonance system or adjacent rooms to a magnetic resonance system.


Note: Sites should ensure that their security systems and other equipment do not interfere with Fenom Pro.

## Chapter 3: Fenom Pro Quick Start Guide

To perform a FeNO test, follow these three simple steps. For full test guidelines and instructions, see Chapters 4 and 5. For setup instructions, see Chapter 9.


**NOTE:** Check to make sure the device is powered on. If the device is on but the display is blank, touch the screen to wake up the device. The device may take one minute to warm up.

1. Select **Start** test on the **Home Screen**. Remove a new single-patient-use mouthpiece from its packaging and attach it to the hand-piece by pressing the mouthpiece towards the top of the hand-piece and twist clockwise to secure. Be careful to not touch the part of the mouthpiece that will have patient contact.

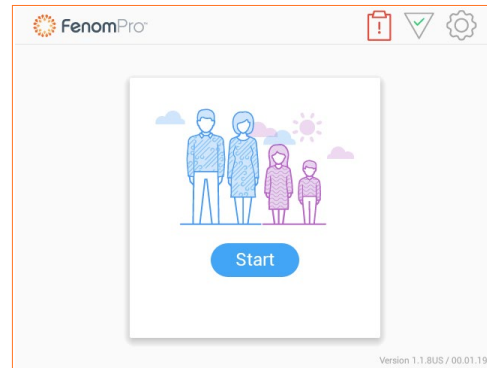


**NOTE:** See Chapter 8: Practice Mode if patient requires a demonstration before taking the test.

2. Touch the **Begin Test** button and instruct the patient to inhale naturally to full capacity, place mouth on the mouthpiece ensuring to keep tight seal so no air escapes and exhale for the full breath test duration at a steady flow.

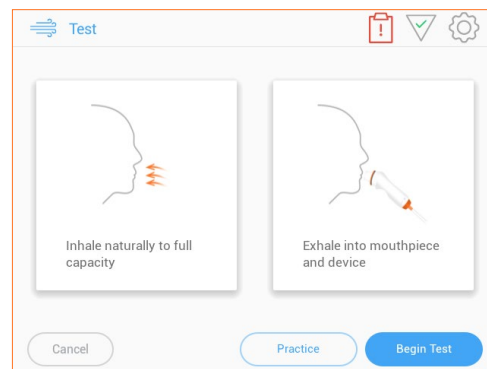


Instruct the patient to keep the **indicator** over the star at the top of the gauge.


**NOTE:** Having the indicator within the green range is also acceptable.

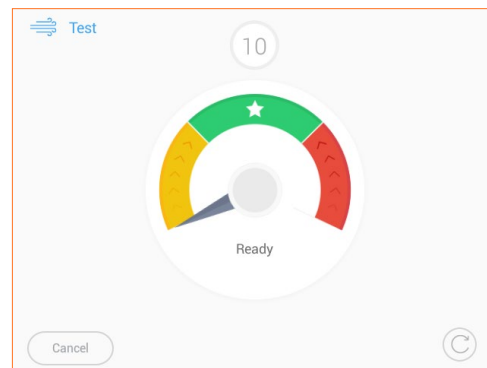


3. The Fenom Pro will display the **Stop Now** screen and play an audible chime once the patient has successfully completed the breath maneuver.

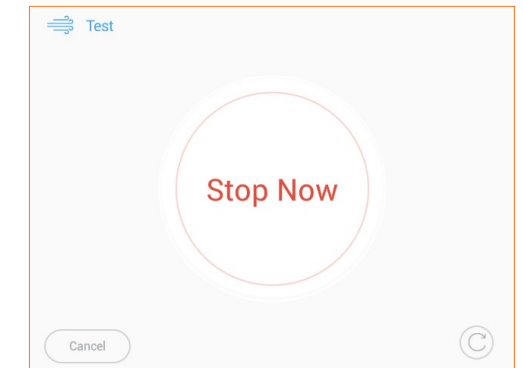

Results will display in 28 seconds.

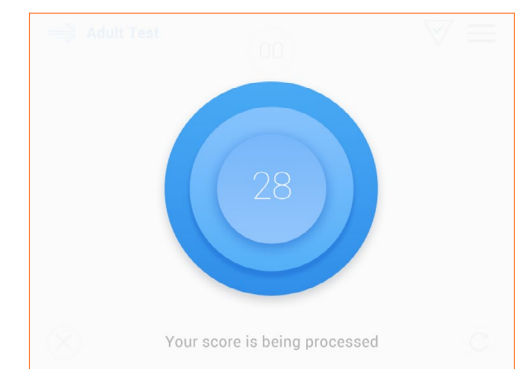

4. Press the **Done** button and properly dispose of the used mouthpiece.

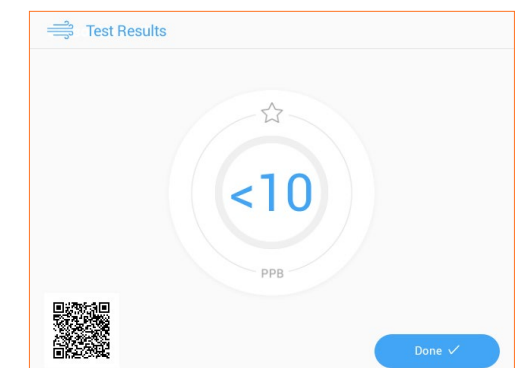

## Chapter 4: FeNO Measurement Preparation


**NOTE:** See *Chapter 2: Safety and Warnings* for list of safety instructions and warnings.

### Wake up Device

1. If the device is powered off, press the Power button to turn it on.

**NOTE:** If the device is powered on but the display is blank, touch the screen to wake up the device.

Allow the device to warm-up for one minute.

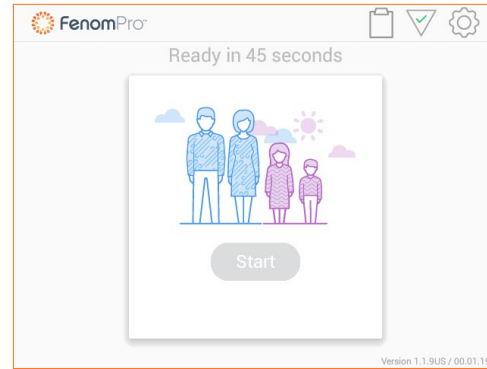

### Pre-test Check

1. Check the battery indicator to ensure the unit has sufficient battery power to perform a FeNO measurement. If the battery indicator is below 25%, plug device into the power supply before using.
2. Check that device is on a flat, stable surface while performing a FeNO measurement.
3. Confirm that the patient meets eligibility requirements:
  - Age 7 and up.
  - Has not consumed food or fluids other than water in the preceding 60 minutes.
  - Has not exercised or smoked in the preceding 60 minutes.
4. When pre-test check is complete, proceed with *Chapter 5: Perform FeNO Measurement*.

## Chapter 5: Perform FeNO Measurement

The FeNO measurement is performed by the patient blowing into a single-patient-use mouthpiece that is attached to a handpiece. The patient must blow into the mouthpiece at a controlled rate, which is monitored through an animated graphic on the touch display. Once a sufficient amount of the patient's breath is captured, the sensor analyzes the breath and reports a FeNO score in parts per billion (ppb).

### Perform a FeNO Test

**NOTE:** Complete the steps in *Chapter 4: FeNO Measurement Preparation* before continuing with the steps below.

1. Lift the handpiece out of the cradle on top of the Fenom Pro.
2. Remove the new, single-patient-use mouthpiece from its packaging without touching the part that will go into the patient's mouth.
3. Attach the mouthpiece to the handpiece by firmly grasping the outer diameter of the mouthpiece and pushing towards the top of the handpiece while twisting clockwise until secure.



4. Hand the handpiece to the patient with the mouthpiece attached.
5. Press the **Start** button on the test selection screen.
6. Provide the patient with a brief overview on how to use the Fenom Pro device.
  - Instruct patient to inhale naturally to full capacity before placing mouth on the mouthpiece.
  - Instruct patient to place mouth on the mouthpiece and exhale for a full **10** seconds at a steady flow.
  - Instruct patient to keep lips sealed around the mouthpiece so no breath escapes from patient's lips.


**NOTE:** See *Chapter 8: Practice Mode* if patient requires a demonstration before taking the test.

7. Touch the **Begin Test** button when the patient understands the instructions and is ready to begin.
8. The visual incentive gauge displays.
9. Instruct the patient to begin exhaling into the mouthpiece whenever ready.
10. Ensure that the patient stops exhaling once the **Stop Now** screen is displayed.
11. If the patient was unsuccessful in performing a breath maneuver, review the reason for the failure. If necessary, patient can attempt a test in Practice Mode (Chapter 8) before repeating.
12. Proceed to instructions in *View Results* section of this Chapter.

## View Results

Upon completion of the FeNO test, the patient's breath is analyzed, and the results are displayed in ppb. It takes approximately 28 seconds for the results to display.

1. View the FeNO result. If the result is less than 10 ppb, "<10" will be displayed. If the result is greater than 200 ppb, ">200 ppb" will be displayed.
2. Touch the **Done** button.

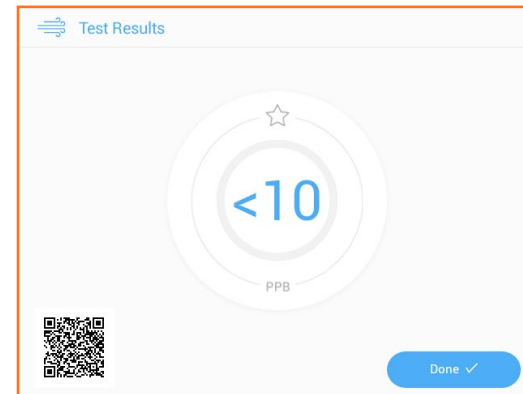

## Remove Mouthpiece

When the patient has completed performing a FeNO measurement:

1. Remove the mouthpiece by firmly grasping around the outer diameter and twist counter-clockwise while pulling away from the handpiece.
2. Properly dispose of the used mouthpiece.
3. Replace the handpiece in its cradle on top of the device.



## Chapter 6: Quality Control

The quality control (QC) measurement is performed on a Fenom Pro device by qualified operators. QC Mode is designed to ensure the instrument is operating within its specifications. The QC procedure must always be activated as a daily QC measurement when the device is used in a clinical setting. The device will display the results of QC testing, users are only required to read Passed, Failed, or Expired result.

To be a valid QC User, a health care professional shall satisfy the following criteria: over 18 years of age, non-smoker, no known airway disease or chronic cold, preferably no allergies or asthma, and stable FeNO values between 10 and 40 ppb.

### Perform a QC Test

1. Touch the QUALITY CONTROL STATUS button from the Settings screen or the title bar at the top of the screen.
2. Touch the SELECT QC USER button.
3. Select the user listed in the SELECT USER FOR TEST window that corresponds to the actual user that will perform a test. If the user is not listed, create a new user by selecting the NEW USER button and following the instructions in the Create QC User section of this chapter.
4. Once a user is selected, touch SELECT TEST button.
5. Follow the instructions in Chapter 5: Perform FeNO Measurement to perform a QC test

**NOTE:** Each day, the system requires a valid QC user and a Negative Control User to perform a QC test. QC tests for the Negative Control User can be performed by any QC user. The Negative Control User only needs to perform the breath maneuver to trigger the measurement. A negative control consisting of a NO free gas sample is generated by the device after the Negative Control User activates the measurement by successfully completing the breath maneuver.

### QC User Status

**NOTE:** Four tests are required for a given QC user in order to establish quality control limits. If a QC User result fails, then an alternative QC User should perform the QC Measurements to determine if the first QC User may have a changing condition or if the system may require service. If the 2nd QC User fails QC please contact customer support.

- **Passed** - latest test is within the expected range for QC user.
- **Failed** - latest test is outside the expected range for QC user.
- **Expired** - 24 hours or greater have elapsed since the last passing test for the QC user.
- **Insufficient Data** - fewer than four tests have been performed for the QC user.

### Create QC User

1. Touch NEW USER button from the SELECT USER FOR TEST window.
2. Enter USERNAME. The USERNAME may not contain spaces or match existing users.
3. Touch SAVE button.



## Chapter 7: Power Off Device

It is OK to leave the device powered-on at all times. Your device will automatically go into sleep mode when not being used. Only power off if you don't intend to use for extended periods of time.

### To power off the device:

1. Hold the **POWER** button down for at least 1 second.
2. Touch **OK** on the confirmation window.

**NOTE:** It is recommended to keep the Fenom Pro device connected to a power supply whenever possible when powered on.

## Chapter 8: Practice Mode

Practice mode is an option for use by a new patient in order to demonstrate the steps for performing a FeNO test. Results are not recorded in this mode.

### IMPORTANT!

Using the Practice mode counts as one exhalation toward the maximum six exhalations per patient per day and toward the maximum three exhalations per mouthpiece

### To access the Practice mode:

1. On the **Home** screen, select the **Start**.
2. Touch the **Practice** button.
3. Attach a new, single-patient-use mouthpiece to the handpiece and review how to use the Fenom Pro device with the patient. (See *Chapter 5: Perform FeNO Measurement* for detailed instructions.)
4. The patient inhales then begins exhaling into the mouthpiece for 10 seconds when ready.
5. The patient stops exhaling once the countdown reaches 0 (zero).
6. If training was successful, **Good job!** displays.
7. If training was unsuccessful, **Try again** displays.
8. Touch **Repeat** arrow to perform training again and go back to Step 4, or touch **Done** if finished to go back to the **Home** screen.

## Chapter 9: Device Setup

### Initial Setup

#### To set up the Fenom Pro device:

1. Remove the device and power cable from the shipping package.

**NOTE:** Retain all packaging for future transportation of the device.

2. Connect the breath tube to the orange port at the bottom of the handpiece. Ensure the breath tube is fully seated against the rear surface, as shown in the figure. The second figure shows an improper breath tube connection.



Correct



Improper

3. Once the breath tube is connected, place the handpiece in the cradle on top of the device.
4. Connect the power cable from the rear panel of the device to an outlet. (See *Table 1* for power connection location.) The AC Power Indicator displays green when the device is plugged in and powered on. **NOTE:** The device should be allowed to charge for at least 4 hours before operating on battery power. The device can operate normally while charging.
5. Press the **Power** button to turn on the device.
6. The Device Setup screen will display after powering on. From this screen, set the following device settings:
  - Select Language
  - Set Time/Date
  - Add Tests


These settings can be accessed and changed at any time.
7. The Fenom Pro is now ready to begin a test.

## Configuration Settings

There are device settings that require configuring based on location and requirements. These settings are accessed through the Settings icon. (See *Table 2: Indicators and Icons.*)


### Device Settings

1. The Settings button provides access to set the Time and Date, select Language, view System Info, select Test Incentive Sound level, and Add Tests.

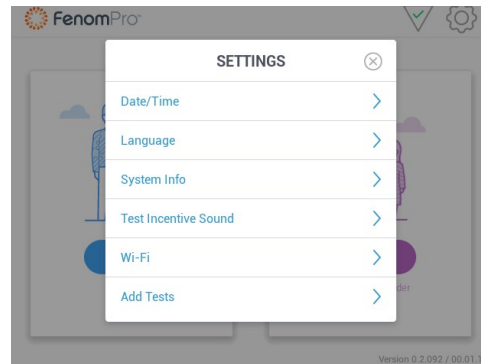


### Time/Date

1. Touch Time/Date on the Settings screen to set the date and time on the device.
2. Use the + and - buttons to set the date and time.
3. Touch the AM/PM button to toggle between values.
4. Touch the Time Zone drop-down list and select the correct time zone.

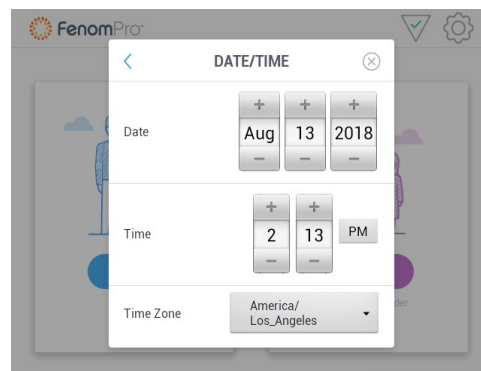


### Language

1. Touch Language on the Settings screen.
2. Select the desired language.
3. A check mark next to the language indicates the selected language.

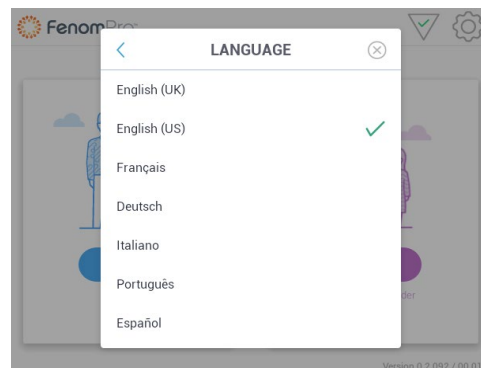


### System Information

1. Touch System Info on the Settings screen to view Device Serial Number, Licensed Tests (number of licensed tests remaining), Service Due Date, Software Version, and Firmware Version.

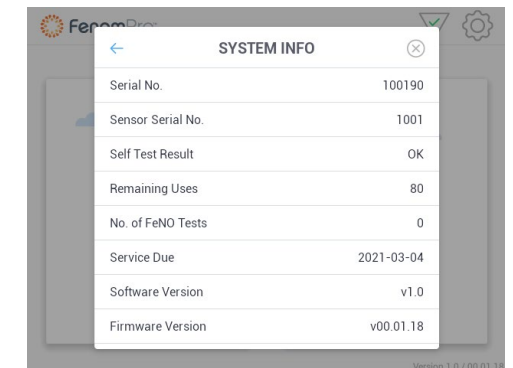


### Test Incentive Sound

1. Touch Test Incentive Sound on the Settings screen to set the desired volume level, High, Low, or Off.
2. A check mark indicates the current selection.

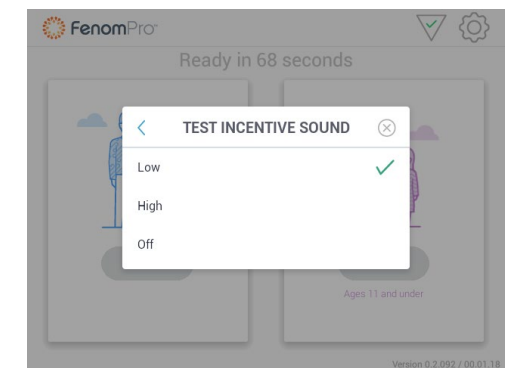


### Add Tests

1. Contact Spirosure support and ask for additional licensed tests for the Fenom Pro.
2. Enter the code provided by the support representative and press Add Tests.

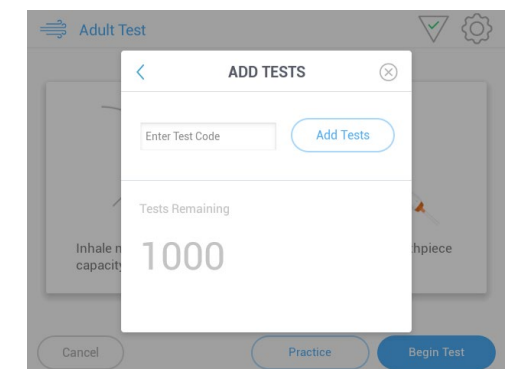

### Quality Control Status

1. Touch Quality Control Status on the Settings screen.
2. See Chapter 6 Quality Control for more instructions..

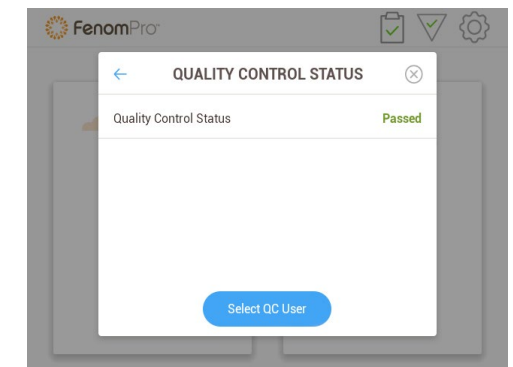

## Chapter 10: General Care

Follow the recommendations below for cleaning and general care of the Fenom Pro and its accessories.

### IMPORTANT!

Never attempt to open or service the Fenom Pro device or component.

### Operating Conditions

Ensure stable operating conditions by avoiding placement of the device in direct sunlight, near sources radiating heat, or ventilation. The device operates under the following conditions:

- Temperature range of 15 to 30°C (59 to 86°F)
- Atmospheric pressure range of 106 to 80 kPa
- Relative humidity range of 20 to 80%, non-condensing

### Cleaning

- Clean the external surfaces of the device with a cloth pre-moistened with 5% bleach solution at the end of each day of use.
- DO NOT use spray detergents.

### Handling

- Take care while handling the device.
- DO NOT drop the device or the handpiece.
- Carry the device by placing fingers in the recessed handle on the back and placing thumb over the top of the device. Support the device from the bottom with other hand.

### Storage

- Clean the device before storing.
- Store the device in its original shipping packaging.
- Store the device in a location free from dust, free from excessive moisture or water splash, and away from excessive heat, cold, or dry conditions.
- DO NOT store the device on tall or unstable surfaces.
- Store mouthpieces in original, unbroken packaging.

### Preventive Inspections

- Ensure the handpiece is not damaged and is in good condition.
- Ensure the tubing from the handpiece to the unit is not damaged and is in good condition.
- Ensure the power cord is not damaged and is in good condition.
- Ensure the touch screen is not damaged and is in good condition.

### Rechargeable Battery

- Use only the power adapter provided by Spirosure to charge the Fenom Pro device.
- Capacity: > 15 tests over 6 hours on a fully charged battery
- Charging time: 4 hours

Lowered capacity: Extended charge times or reduced operation indicates the battery should be replaced. Contact Spirosure support for assistance.

### Maintenance

- Periodic service is required. Check system information (page 19) for service due date.
- Contact Spirosure support to schedule service.

### Disposal of Used/Expired Equipment and Consumables

- Expired devices should be recycled according to the local program for electronic equipment.
- Used or expired mouthpieces should be recycled according to the local program.



### Limited Warranty

Spirosure, Inc. warrants the Fenom Pro to be free of defects in materials and workmanship for a period of 18 months from the date of shipment. Spirosure's sole obligation under this warranty is limited to repairing or replacing, at its choice, any item covered under this warranty when such an item is returned intact and prepaid, to Spirosure or the local representative.

The product warranty is automatically invalidated if the products are repaired, altered or otherwise tampered with by unauthorized personnel, or have been subject to misuse, neglect or accident.

The product warranty does not cover product failure or damage resulting from use with non-approved accessories. Spirosure takes no responsibility for health and safety problems or other problems caused by the use of accessories not approved by Spirosure.

# Chapter 11: Troubleshooting

The Fenom Pro device, subcomponents, and accessories are not field serviceable.

## Support

Please contact Spirosure support if the Fenom Pro presents any problems that cannot be solved with the actions stated in this manual.


## Add Tests

The Fenom Pro requires licensed tests in order to perform FeNO measurements. When the number of licensed tests approaches zero, the Licensed Tests Status Button will turn red.

### To order additional tests, follow three simple steps.

1. Contact Spirosure support representative and ask for additional licensed tests for the Fenom Pro.
2. Press the Licensed Tests Status button or navigate to the Add Tests screen in the Settings menu.
3. Enter the code provided by the support representative and press Add Tests.

**NOTE:** Please contact Spirosure support if the Fenom Pro device is approaching the service date or zero remaining licensed tests.

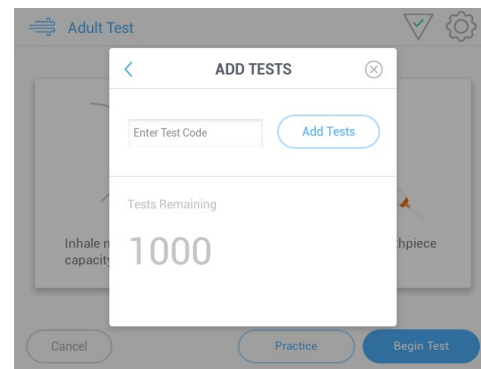

## Error Codes

In the event the device displays an error message, use Table 4 to look up the error code and perform the suggested actions to resolve the issue.

Table 4: Error Codes

Error Code	Error Situation/Error Message	Actions
	Power up self-test fails.	Cycle the power on the device. After two consecutive self-test failures, contact Technical Support.
10-011	Please contact customer support to continue using this device. FeNO tests are disabled.	Warning that device has reached the maximum number of uses or the expiration date; contact Technical Support.
10-012	Please contact customer support to continue using this device. FeNO tests are disabled.	Test not allowed due to device having reached the maximum number of uses or the expiration date; contact Technical Support.
10-013	Device time is incorrect. Please set device time in Settings.	See <i>Chapter 8: Device Settings – Configuration Setting</i> section for instructions on how to set the time.
10-014	Please contact customer support to continue using this device. FeNO tests are disabled.	Contact Technical Support. Test not allowed due to device having reached the maximum warm hours.
10-019	The test was stopped because breath airflow fell below the minimum threshold. Please try again.	Give the patient a moment to rest, restate the proper breath maneuver pointing out the <b>green</b> target zone, and then try test again. If problem persists, try performing a shortened test. (See <i>Chapter 10: Shortened Test Mode (For Research Only)</i> .)
10-020	The test was stopped because breath airflow exceeded the maximum threshold. Please try again.	Restate the proper breath maneuver pointing out the <b>green</b> target zone, and then try the test again.
10-021	The test was stopped because breath airflow was out of the desired range for too long. Please try again.	Give the patient a moment to rest, restate the proper breath maneuver pointing out the <b>green</b> target zone, and then try the test again. If problem persists, try performing a shortened test. (See <i>Chapter 10: Shortened Test Mode (For Research Only)</i> .)
10-025	The test was stopped because breath airflow started too early or continued after the breath maneuver. Please try again.	Restate the proper breath maneuver pointing out the <b>green</b> target zone, and then try the test again.
10-043	Battery level is very low. Please immediately plug in the device or turn it off.	Plug device into power supply before using.
10-064	Battery level has fallen below what is required to do a test. Please plug in the device.	Plug device into power supply before using.
10-065	Warning that the device sensor has reached the expiration date. FeNO tests are disabled.	Please contact customer support to continue using this device.

Error Code	Error Situation/Error Message	Actions
10-066	Warning that the device date is incorrect	Set the device date in Settings.
20-001	There was a failure on the device communicating with the sensors. If problem persists contact support. Please restart the device and try again.	Cycle the power on the device and try again. If problem persists, contact Technical Support.
20-002	There was a failure on the device. Please restart the device and try again.	Cycle the power on the device and try again. If problem persists, contact Technical Support.
40-028	There was an error calculating reading. Please wait a few minutes and try again.	Wait a few minutes and try again. If problem persists, contact Technical Support.
40-065	The test was stopped because the pump airflow was outside the allowed threshold. Please try again.	Wait one minute and try again. If problem persists, contact Technical Support.
40-066	The test was stopped because flow variability was outside the allowed threshold. Please try again.	Wait one minute and try again. If problem persists, contact Technical Support.
40-067	The test was stopped because peak readings were beyond maximum. Please try again.	Wait one minute and try again. If problem persists, contact Technical Support.
40-068	The test was stopped because baseline readings were below minimum. Please try again.	Wait one minute and try again. If problem persists, contact Technical Support.
45-044	A failure has occurred on the device hardware. Please contact support. Issue: Unknown failure code.	Contact Technical Support.
45-045	A failure has occurred on the device hardware. Please contact support. Issue: Memory CRC Error.	Contact Technical Support.
45-048	A failure has occurred on the device hardware. Please contact customer support. Issue: Battery Communication Failure.	Contact Technical Support.
45-050	A failure has occurred on the device hardware. Please contact support. Issue: Calibration EEPROM Communication Failure.	Contact Technical Support.

Error Code	Error Situation/Error Message	Actions
45-051	A failure has occurred on the device hardware. Please contact support. Issue: Calibration EEPROM CRC Failure.	Contact Technical Support.
45-052	A failure has occurred on the device hardware. Please contact customer support. Issue: Android Communication Timeout.	Contact Technical Support.
45-054	A failure has occurred on the device hardware. Please contact customer support. Issue Board Temperature High.	Contact Technical Support.
45-063	A failure has occurred on the device hardware. Please contact support. Issue: Battery Charger Failure.	Contact Technical Support.
45-069	A hardware failure has occurred. Please wait while the issue is being resolved. If the issue persists please restart the device.	Wait a few minutes for issue to resolve. If problem persists, cycle the power on the device and try again. If problem persists, contact Technical Support.
90-040	A failure has occurred on the device license. FeNO tests are disabled.	Please contact customer support to continue using this device.

## Chapter 12: Technical Data

<b>Dimensions and Weight</b>	Height: 145 mm Width: 230 mm Depth: 140 mm Weight (including handpiece): 2.4 kg
<b>Electrical Data</b>	Device power consumption: < 20 VA Power supply mains voltage: 100–240 V ~50–60 Hz
<b>Exhaled NO Performance</b>	The Fenom Pro is verified to fulfill performance herein under temperature range of 15–30°C (59–86°F), relative humidity of 20–80%, and pressure range of 106–80 kPa.
<b>Linearity</b>	Slope 1.00 ± 0.05 Squared correlation coefficient, r2 ≥ 0.998
<b>Precision</b>	NO concentrations ≤ 50 ppb: 5 ppb NO concentrations > 50 ppb: 10% of the concentration
<b>Accuracy</b>	NO concentrations ≤ 50 ppb: ± 5 ppb NO concentrations > 50 ppb: ± 10% of the concentration
<b>Limit of Detection</b>	10 ppb
<b>Measurement Range</b>	10–200 ppb
<b>Exhalation Parameters</b>	Exhalation time: 10 seconds Exhalation pressure is between 15–20 cm (6–8 in.) water Exhalation flow rate is 45–55 mL/s; warning sounds played outside of this range

## Chapter 13: Device Performance

### Linearity

Two Fenom Pro devices were tested as part of the linearity study. Nitric oxide (NO) was mixed in a balance gas of air to obtain 1 NO concentration levels between <10 and >300 ppb. The data showed that the test is linear across this range. The slope of the fitted line was 1.0 ± 0.05 and the r2 coefficient was > 0.998. The results are listed in the table below.

	Slope	Intercept	R2
Combined Device	1.05	-2.16	1.00

### Analytical Precision

Analytical precision was evaluated using samples across 10, 25, 75 and 200 ppb concentrations. Two replicate determinations of each concentration were made eight times a day for five days across five different devices.

Repeatability is an estimated of variation withing one test run in one day. Within-device precision is an estimate of variation between test runs and days. The results for 10 and 25 ppb are expressed as absolute values in ppb. The results for the 75 and 200 ppb levels are expressed as percentage of the measured NO concentrations. Both standard deviation estimates meet the precision claim at all four concentration levels. The results are listed in the table below.

Device	Repeatability				Within-Device Precision			
	10 ppb	25 ppb	75 ppb	200 ppb	10 ppb	25 ppb	75 ppb	200 ppb
100183	1.1	1.5	2.4%	1.4%	2.0	3.2	3.7%	1.3%
100184	0.8	1.3	2.1%	1.4%	1.0	1.9	3.5%	1.5%
100186	0.8	1.2	1.8%	1.1%	1.5	2.0	3.8%	1.3%
100187	1.2	1.7	5.0%	4.5%	1.5	2.5	5.5%	4.5%
100193	0.8	1.1	1.9%	1.7%	1.4	1.5	2.4%	1.8%

### Accuracy

Four Fenom Pro devices were tested over five concentrations between 10 and 200 ppb NO with five replicates each across five different environmental conditions. Fenom Pro results were compared to results obtained from an independent device calibrated against NIST values. The data showed that mean Fenom Pro results were within 5 ppb for samples less than 50 ppb, and percent errors were less than 10% for samples at 50 ppb and greater. The results are listed in the table below. The accuracy for all test cases within the technical specification.

Device	10 ppb	25 ppb	75 ppb	100 ppb	200 ppb
100179	-1.8	-2.8	-3.8%	-1.4%	0.4%
100180	-2.0	-2.2	-0.4%	0.4%	3.4%
100183	-2.0	-3.4	-1.8%	0.0%	2.8%
100184	-1.8	-1.6	-0.8%	-0.2%	1.6%

### Limit of Detection

Two Fenom Pro devices were tested using NO gas with air balance at 3 and 5 ppb. Five replicate determinations of each concentration were made. The mean and confidence limits at 3 ppb and 5 ppb support that the detection limit is below 10 ppb, supporting the claim. Results are listed in the table below.

Device		3 ppb	5 ppb
100175	Mean	3.9 ppb	5.1 ppb
	Lower 95% CI	3.5 ppb	4.7 ppb
	Upper 95% CI	4.3 ppb	5.5 ppb
100177	Mean	3.9 ppb	5.2 ppb
	Lower 95% CI	3.3 ppb	4.6 ppb
	Upper 95% CI	4.5 ppb	5.8 ppb

## Chapter 14: Reference

### Symbol Explanation

	WEEE Directive 2102/19/EU		Manufacturer
	Keep out of rain & damp conditions		Do not reuse
	Prescription Only		In vitro diagnostic device
	MR Unsafe – Fenom Pro is not rated for use near magnetic resonance		Use by YYYY-MM-DD (expiry)
	Caution, consult accompanying documents		Catalog part number
	Type BF applied part complying with IEC 60601-1		Lot number
	Non-Sterile		Serial number
	Consult instructions for use		Quantity
	Operating humidity range		Operating temperature range
	Intertek ETL Listed, Canada & USA		

## Chapter 15: Parts and Accessories

### Warning!

Any accessory not recommended by Spirosure, Inc. may result in loss of performance, damage to your Fenom PRO, or injury. The product warranty does not cover product failure or damage resulting from use with non-approved accessories. Spirosure, Inc. takes no responsibility for health and safety problems or other problems caused by the use of accessories not approved by Spirosure.

### Parts

- Fenom Pro Model No. 900-0001
- Packaged Fenom Pro P/N: 900-0005
- Fenom Pro Power Cable P/N: 415-0003
- Fenom Pro Power Supply P/N: 197-0001
- Fenom Pro IFU Reference P/N: FRM-3301

### Accessories

- Fenom Pro Single Patient Use Mouthpiece\* (20 Count) P/N: 900-0019
- \* Disposable mouthpiece to be changed for every patient.

To order accessories, contact Spirosure support. If you are unsure of your representative's contact information, visit [www.fenomasthma.com](http://www.fenomasthma.com).

## Bibliography

- [1] "ATS/ERS Recommendations for Standardized Procedures for the Online and Offline Measurement Exhaled Lower Respiratory Nitric Oxide and Nasal Nitric Oxide, 2005." Am. J. Respir. Crit. Care Med., 2005; vol. 171, pp. 912–930.
- [2] A. H. Henriksen, M. Sue-Chu, T. L. Holmen, A. Langhammer, and L. Bjermer, "Exhaled and nasal NO levels in allergic rhinitis: relation to sensitization, pollen season and bronchial hyperresponsiveness." Eur. Respir. J., 1999; vol. 13, pp. 301–306.
- [3] G. Rolla, L. Brussino, M. T. Bertero, P. Colagrande, M. Converso, C. Bucca, S. Polizzi, and F. Caligaris-Cappio, "Increased nitric oxide in exhaled air of patients with systemic lupus erythematosus." J. Rheumatol., 1997; vol. 24, pp. 1066–1071.
- [4] C. Soderman, A. Leone, V. Furst, and M. G. Persson, "Endogenous nitric oxide in exhaled air from patients with liver cirrhosis." Scand. J. Gastroenterol., 1997; vol. 32, pp. 591–597.
- [5] D. Cordeiro, A. Rudolphus, E. Snoey, and G. J. Braunstahl. "Utility of nitric oxide for the diagnosis of asthma in an allergy clinic population." Allergy Asthma Proc., 2011; vol. 32, pp. 119–126.
- [6] K. Ansarin, J. M. Chatkin, I. M. Ferreira, C. A. Gutierrez, N. Zamel, and K. R. Chapman, "Exhaled nitric oxide in chronic obstructive pulmonary disease: relationship to pulmonary function." Eur. Respir. J., 2001; vol. 17, pp. 934–938.
- [7] K. Alving and A. Malinowski, "Basic aspects of exhaled nitric oxide." Exhaled Biomarkers, I. Hovarth and J. C. de Jongste, Eds. European Respiratory Society, 2010, pp. 1–31.
- [8] A. C. Olin, A. Aldenbratt, A. Ekman, G. Ljungkvist, L. Jungersten, K. Alving, and K. Toren, "Increased nitric oxide in exhaled air after intake of a nitrate-rich meal." Respir. Med., 2001; vol. 95, pp. 153–158.
- [9] R. J. Meijer, H. A. Kerstjens, D. S. Postma, G. H. Koeter, T. W. van der Mark, "Exhaled nitric oxide concentration is influenced by alcohol containing disinfectants." Eur. Respir. J., 1996; vol. 9, p. 1111.
- [10] R. A. Dweik, P. B. Boggs, S. C. Erzurum, C. G. Irvin, M. W. Leigh, J. O. Lundberg, A. C. Olin, A. L. Plummer, D. R. Taylor, "An Official ATS Clinical Practice Guideline: Interpretation of Exhaled Nitric Oxide Levels (FENO) for Clinical Applications." Am. J. Respir., 2011; vol. 184, pp. 602–615.



**Spirosure, Inc.**  
 7020 Koll Center Pkwy, Suite 110  
 Pleasanton, CA 94566-3107 USA  
 (888) 609-4839  
 support@spirosure.com

